


A New Insufflator & Smoke Evacuation System

EVA15 is a new insufflator and smoke evacuation system for safe and smoke-free operating rooms. The compact and pedal-activated EVA15 smoke evacuation system safely filters the hazardous contaminants of surgical smoke before they become airborne in the OR. EVA15 compensates for leaks and suctioning, while delivering a constant insufflation pressure down as far as 7mmHg.



EVA15

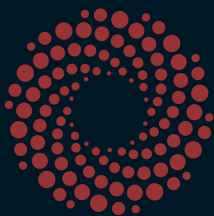

Surgical smoke can no longer be accepted as a normal part of the OR environment.

Product Features

- Works with any standard laparoscopic trocar
- Continuous 7-15mmHg pressure insufflator compensates for leaks
- Pedal activated smoke evacuation
- Rapid end of case desufflation
- ULPA activated carbon evacuation filter (99.9997% efficiency at 0.04 micron)
- Meets AORN guidelines for surgical smoke safety
- Connects to OR scavenging system
- Simple intuitive user interface

Description	Part Number
EVA15 Insufflator and Smoke Evacuation System	EVA15
3 line tubeset for laparoscopic insufflation	LAP-300

